

CPAP Masks
CPAP Machines
AUTO Machines
BPAP Machines



Sleep Well, Live Better



Certified CPAP Supplier



Sleep Health FOUNDATION

CODE OF PRACTICE

FOR SUPPLIERS OF CPAP SLEEP THERAPIES

CPAP SALES[®]

Making sleep apnea healthcare more affordable

Customer Support

Phone: 1300 04 CPAP

Ph: 1300 04 2727
Fax: (02) 4744 2429
Email: support@cpapsales.com.au
Website: www.shop.cpapsales.com.au

PO Box 118
Riverstone
NSW 2765
Australia

www.shop.cpapsales.com.au

CPAP SALES®

Making Sleep Apnea Healthcare More Affordable



RESmart CPAP and Auto Machines with Heated Humidifier



| Models | Mask | Function | Price \$ |
|-----------|----------|-----------------------|------------------|
| CP10H | No | CPAP (Fixed Pressure) | Refer to website |
| CP10BONUS | Included | CPAP (Fixed Pressure) | Refer to website |
| AP10H | No | CPAP/AUTO/TITRATE | Refer to website |
| AP10BONUS | Included | CPAP/AUTO/TITRATE | Refer to website |

Luna GII CPAP and Auto Machines with Heated Humidifier



| GII Models | Mask | Function | Price \$ |
|-------------|------|-----------------------|------------------|
| CPAP2GH/2.4 | No | CPAP (Fixed Pressure) | Refer to website |
| APAP2GH/2.4 | No | CPAP/AUTO/TITRATE | Refer to website |
| APAP2GH/3.5 | No | CPAP/AUTO/TITRATE | Refer to website |

BPAP – Bi-level Machines with Heated Humidifier



| Models | Pressure | Function | Auto | Price \$ |
|----------|-----------|-----------------|------|------------------|
| BPAP 25 | 4~25cmH2O | S, CPAP | No | Refer to website |
| BPAP 25A | 4~25cmH2O | S, CPAP | Yes | Refer to website |
| BPAP 25T | 4~25cmH2O | S, CPAP, T, S/T | No | Refer to website |
| BPAP 30T | 4~30cmH2O | S, CPAP, T, S/T | No | Refer to website |

CPAP Masks - Nasal, Nasal Pillow and Full Face



| CPAP Masks | Sizes | Price \$ |
|--------------|-----------------------------|------------------|
| Nasal Pillow | S / M L pillows included | Refer to website |
| Micro Nasal | S or M or L sizes available | Refer to website |
| Full Face | S or M or L sizes available | Refer to website |

For the latest prices, please visit www.shop.cpapsales.com.au

CONTACT CPAP SALES PTY LTD
 PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au

CPAP SALES®

Making Sleep Apnea Healthcare More Affordable



Transcend II Travel CPAP, Auto Machines and Accessories



| Transcend CPAP Options | Mask | Function | Price \$ |
|-------------------------------|------|----------|------------------|
| CPAP Starter Kit | No | CPAP | Refer to website |
| CPAP Starter Kit | Yes | CPAP | Refer to website |
| CPAP Starter Kit & P4 Battery | No | CPAP | Refer to website |
| CPAP Starter Kit & P4 Battery | Yes | CPAP | Refer to website |
| CPAP Starter Kit & P8 Battery | No | CPAP | Refer to website |
| CPAP Starter Kit & P8 Battery | Yes | CPAP | Refer to website |



| Transcend Auto Options | Mask | Function | Price \$ |
|-------------------------------|------|-----------|------------------|
| Auto Starter Kit | No | CPAP/Auto | Refer to website |
| Auto Starter Kit | Yes | CPAP/Auto | Refer to website |
| Auto Starter Kit & P4 Battery | No | CPAP/Auto | Refer to website |
| Auto Starter Kit & P4 Battery | Yes | CPAP/Auto | Refer to website |
| Auto Starter Kit & P8 Battery | No | CPAP/Auto | Refer to website |
| Auto Starter Kit & P8 Battery | Yes | CPAP/Auto | Refer to website |



| Transcend Accessory Options | Price \$ |
|-----------------------------|------------------|
| P4 – 1 Night battery | Refer to website |
| P8 – 2 Night battery | Refer to website |
| Mobile Power Adaptor | Refer to website |
| Docking Base Station | Refer to website |
| Heated Humidifier | Refer to website |

CONTACT CPAP SALES PTY LTD
 PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au

Comfort and Cleaning, Parts and Accessories



| Comfort and Cleaning | | Model # | Price \$ |
|-------------------------------|--|-------------|------------------|
| CPAP Pillow | | CPAP-Pillow | Refer to website |
| Antibacterial CPAP wipes | | CMW50 | Refer to website |
| CPAP Hose Cover | | CHC | Refer to website |
| CPAP Heated Hose | | HH10G | Refer to website |
| CPAP Hose | | SH10G | Refer to website |
| Seatec Chin Strap S/M/L | | CSS/SML | Refer to website |
| Luna Accessories | | Model # | Price \$ |
| Air filter for Luna 5 pack | | LAF5 | Refer to website |
| Easy Clean Water Tank | | LWT | Refer to website |
| Spo2 Kit for Luna CPAP / Auto | | SP02KIT | Refer to website |
| Carry Bag For Luna | | CB10G | Refer to website |
| RESmart Accessories | | Model # | Price \$ |
| Air filter for RESmart 5 pack | | AF05G | Refer to website |
| Air filter cover | | FC10A | Refer to website |
| Easy clean water tank | | YT2A | Refer to website |
| Carry Bag for RESmart | | CB10G | Refer to website |
| Software and Data cable | | SW10GLDC | Refer to website |
| Humidifier with water tank | | HU10A | Refer to website |



CPAP Machine with Heated Humidifier

The BMC RESmart fixed pressure CPAP machine has been designed to provide quality, affordable treatment of Obstructive Sleep Apnea (OSA) in home and clinical use.



Features include...

- CPAP operation (4~20cmH20)
- Exhalation pressure relief (0~3)
- 0 ~ 60 minute ramp time
- Removable heated humidifier
- Powerful on-board memory with iCode® (No SD Card required)
- Easy review of the last night's sleep data
- Automatic airflow start when breathing detected
- Display with backlight and illuminated user buttons
- Convenient carry handle
- Carry Bag, Instructions and spare air filter included

**Refer to
website for
latest
pricing**



BMC MASK OF YOUR CHOICE (normally \$169~\$199)

UPGRADE TO 2 YEAR WARRANTY (normally \$100)

SLEEP THERAPY SOFTWARE (normally \$110)

1.8M DATA CABLE (normally \$39)

UNLIMITED PHONE SUPPORT



quality, affordable, healthcare

ARTG# 218218

CPAP SALES PTY LTD

PO Box 118 Riverstone, Sydney, NSW 2765

Website: www.shop.cpapsales.com.au

Email: sales@cpapsales.com.au

Phone: 0427 299 083 or 1300 04 2727

Fax: (02) 4744 2429

CONTACT CPAP SALES PTY LTD

PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au



AUTO Machine with Heated Humidifier

The BMC RESmart Auto (APAP) Machine includes advanced sensor technology to track your respiratory condition and automatically adjusts to the correct pressure.



Features include...

- Flexible tri-mode operation (CPAP, APAP and Titration)
- Exhalation pressure relief (0~3)
- 0 ~ 60 minute ramp time
- Removable heated humidifier
- Powerful on-board memory with iCode® (No SD Card required)
- Easy review of the last night's sleep data
- Automatic airflow start when breathing detected
- Automatic mask leak detection and compensation
- Display with backlight and illuminated user buttons
- Convenient carry handle
- Carry Bag, Instructions and spare air filter included

**Refer to
website for
latest
pricing**



BMC MASK OF YOUR CHOICE (normally \$169~\$199)

UPGRADE TO 2 YEAR WARRANTY (normally \$100)

SLEEP THERAPY SOFTWARE (normally \$110)

1.8M DATA CABLE (normally \$39)

UNLIMITED PHONE SUPPORT



quality, affordable, healthcare

ARTG# 218218

CPAP SALES PTY LTD

PO Box 118 Riverstone, Sydney, NSW 2765

Website: www.shop.cpapsales.com.au

Email: sales@cpapsales.com.au

Phone: 0427 299 083 or 1300 04 2727

Fax: (02) 4744 2429

CONTACT CPAP SALES PTY LTD

PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au



BPAP (Bi-Level) Machine with Heated Humidifier

This intelligent BPAP machine is a Bi-Level device designed for patients needing different pressure levels for inhalation and exhalation.



Typical uses would be for the clinical therapy of Obstructive Sleep Apnea (OSA), Central Sleep Apnea (CSA) and/or Chronic Obstructive Pulmonary Disease (COPD). The BMC BPAP Series provides Non Invasive Ventilation (NIV) for patients requiring;

- Higher Treatment Pressures (up to 30cmH20)
- Lower Exhalation Pressures

BMC BPAP Machines with Heated Humidifier –

| Model | Pressure | Function | Auto | RRP |
|----------|-----------|-----------------|------|------------------|
| BPAP 25 | 4~25cmH20 | S, CPAP | No | Refer to website |
| BPAP 25A | 4~25cmH20 | S, CPAP | Yes | Refer to website |
| BPAP 25T | 4~25cmH20 | S, CPAP, T, S/T | No | Refer to website |
| BPAP 30T | 4~30cmH20 | S, CPAP, T, S/T | No | Refer to website |

Please consult your respiratory physician for advice on which model is most suitable for you.

Bi-level (BPAP) machines need to be set up strictly in accordance with the settings prescribed by a Respiratory Physician. Regular monitoring and supervision of patients using Bi-Level therapy is required.



CPAP SALES PTY LTD
 PO Box 118 Riverstone, Sydney, NSW 2765
 Website: www.shop.cpapsales.com.au
 Email: sales@cpapsales.com.au
 Phone: 0427 299 083 or 1300 04 2727
 Fax: (02) 4744 2429



Luna GII Auto and CPAP Machines with Heated Humidifier

Introducing the New Generation BMC Luna GII Auto and CPAP Machine.



Large LCD Colour display packed full of information.



New iCode®II for our easiest ever remote reporting of patient compliance data.

Advanced sensor technology to track your respiratory condition and automatically adjusts to the correct pressure.

New Eco Smart dual chamber for more efficient water heating. Opens for cleaning and dishwasher safe.



Features include...

- Flexible tri-mode operation (CPAP, APAP and Titration)
- Exhalation pressure relief (0~3)
- 0 ~ 60 minute ramp time
- Removable Eco Smart heated humidifier
- SD card and iCode® II with free smartphone apps.
- Easy review of the last night's sleep data
- Remote access to sleep therapy data
- SpO2 Pulse Oximeter (optional kit via USB port)



quality, affordable, healthcare
ARTG# 218218

| GII Models | Colour Display | Function | Humidifier | Size (LxWxH) | Weight |
|-------------|----------------|-------------------|------------|--------------|--------|
| CPAP2GH/2.4 | 2.4 inch LCD | CPAP | Yes | 290x196x134 | 2.2kg |
| APAP2GH/2.4 | 2.4 inch LCD | CPAP/AUTO/TITRATE | Yes | 290x196x134 | 2.2kg |
| APAP2GH/3.5 | 3.5 inch LCD | CPAP/AUTO/TITRATE | Yes | 290x196x134 | 2.2kg |

Power Requirements: 100~240VAC 50/60Hz 2.0A (Max) Sound Level: <30dB at 10cmH2O

CONTACT CPAP SALES PTY LTD
PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au



iCode®

iCode® provides remote access to sleep compliance data from any BMC CPAP Machine via a phone call.



Features include...

- iCode® is FREE
 - This feature is enabled on all BMC models
 - iCode® is web based
 - No software installation or direct cable connections are necessary
 - iCode® can be accessed anywhere at anytime
- Sleep clinics, physicians or patients can remotely monitor sleep compliance data.

How does it work?

Sleep compliance data is encoded into a 16 character alpha-numeric string which is displayed on the screen of all BMC machines.

There is one iCode® for each available sleep period.

Enter the serial number of the device plus any or all available iCode® data strings into the free webpage at www.bmc-icode.com

Enter any additional information required to be included in the report such as service provider, patient and physician. Press the Report button to generate a sleep compliance report for printing or emailing.

BMC iCode® – FREE web access for remote access of sleep data



| |
|--|
| <p>CPAP SALES PTY LTD PO Box 118 Riverstone, Sydney, NSW 2765 Website: www.shop.cpapsales.com.au Email: sales@cpapsales.com.au Phone: 0427 299 083 or 1300 04 2727 Fax: (02) 4744 2429</p> |
|--|



iCode® II

iCode® II is the next generation in the remote access and reporting of sleep compliance data.



New iCode® II for our easiest ever remote reporting of patient compliance data.

Features include...

- iCode® II is FREE
 - This feature is enabled on all BMC 2G models
 - iCode® II is web based
 - No software installation or direct cable connections are necessary
 - iCode® can be accessed anywhere at anytime on screen or via the Free Smartphone APP
- Sleep clinics, physicians or patients can remotely monitor sleep compliance data*

How does it work?

Sleep compliance data is encoded into a 16 character alpha-numeric string which is displayed on the screen of all BMC GII machines.

There is one iCode® II for each available sleep period.

Enter the Patient name on the FREE Smartphone APP.

Scan the iCode® II code using the Free Smartphone APP.

View the compliance data from the APP or press the Report button to email the report to your Sleep Clinic, Physician or Health Care Provider.



BMC iCode®II – FREE remote access of sleep data



CPAP SALES PTY LTD

PO Box 118 Riverstone, Sydney, NSW 2765

Website: www.shop.cpapsales.com.au

Email: sales@cpapsales.com.au

Phone: 0427 299 083 or 1300 04 2727

Fax: (02) 4744 2429

CONTACT CPAP SALES PTY LTD


PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au



Micro Nasal CPAP Mask

Dual-wall cushion
soft silicone material
with ergonomic pad
shape provides more
comfort and better seal.

Quickfit clips
easy to fit even in the
dark.

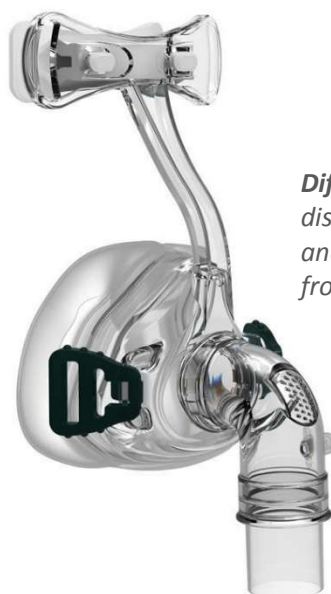

Forehead frame support
Improves comfort and seal.

Diffusing vent
disperses air gently
and quietly away
from bed partner.

360° elbow rotation
positions the tubing
for maximum comfort.

Nasal Pillow CPAP Mask



Adjustable angle of pillows for
best fit.



360° rotating elbow with gentle
and quiet exhaust.



Dual direction
adjustable headgear.



Packed with
S/M/L pillows.



One mask fits most people with
three pillow sizes included.

Super lightweight design provides minimal
face contact for maximum comfort.



quality, affordable, healthcare

ARTG# 219618

All BMC Masks – RRP under \$200

CONTACT CPAP SALES PTY LTD

PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au



Full Face CPAP Mask

More suitable for higher pressures and for 'mouth breathers'



Forehead frame support with a relocatable 'T' bar which can be fitted into three different positions to accommodate a wider range of comfort options.

Bi Lateral diffusing vents disperse air gently and quietly away from bed partner whilst effectively reducing dead space.

Dual-Wall cushion with soft silicone material and an ergonomic pad shape. The surrounding thin pad provides more comfort and a better seal.

Anti-Asphyxia valve provides safety in case airflow is blocked or stops.

360° elbow rotation positions the tubing for maximum comfort.

Mask clips enable the mask to be fitted and removed without the need for continuous re-tensioning of the headgear.

All BMC Masks – RRP under \$200

BMC[®] CPAP
quality, affordable, healthcare
ARTG# 220243

CPAP SALES PTY LTD

PO Box 118 Riverstone, Sydney, NSW 2765

Website: www.shop.cpapsales.com.au

Email: sales@cpapsales.com.au

Phone: 0427 299 083 or 1300 04 2727

Fax: (02) 4744 2429

CONTACT CPAP SALES PTY LTD

PO Box 118 Riverstone, NSW 2765. Ph. 1300 04 2727 www.shop.cpapsales.com.au

CPAP SALES®

Making Sleep Apnea Healthcare More Affordable

BMC® CPAP

- Machine
- Heated humidifier
- Carry bag
- AC power cord
- CPAP hose
- Spare filter
- User manual
- Quick start guide



AUTO

BMC MASK OF YOUR CHOICE (normally \$169~\$199)

UPGRADE TO 2 YEAR WARRANTY (normally \$100)

SLEEP THERAPY SOFTWARE (normally \$110)

1.8M DATA CABLE (normally \$39)

UNLIMITED PHONE SUPPORT

**Fast & Free Delivery
Anywhere in Australia**



PLUS



TRANSCEND™
Sleep Apnea Therapy SYSTEM

CPAP & AUTO

INCLUDES BMC MASK

This material is provided for informational and promotional purposes only and must not be considered as a substitute for professional medical advice, diagnosis or treatment. Any decisions made as a result of viewing this document, or visiting our website should be acted or relied upon only after consultation with your doctor and/or sleep clinic.